# SOEHNLE

## Shape Sense Connect 200

|             |                            |    |
|-------------|----------------------------|----|
| <b>D</b>    | BEDIENUNGSANLEITUNG        | 11 |
| <b>GB</b>   | USER MANUAL                | 14 |
| <b>NL</b>   | GEBRUIKSAANWIJZING         | 17 |
| <b>I</b>    | ISTRUZIONI PER L'USO       | 20 |
| <b>F</b>    | MODE D'EMPLOI              | 23 |
| <b>E</b>    | INSTRUCCIONES DE MANEJO    | 27 |
| <b>P</b>    | MANUAL DE INSTRUÇÕES       | 30 |
| <b>DK</b>   | BRUGSANVISNING             | 33 |
| <b>S</b>    | BRUKSANVISNING             | 36 |
| <b>NO</b>   | BRUKSANVISNING             | 39 |
| <b>FI</b>   | KÄYTTÖOHJE                 | 42 |
| <b>CZ</b>   | NÁVOD K POUŽITÍ            | 45 |
| <b>SK</b>   | NÁVOD NA OBSLUHU           | 48 |
| <b>BG</b>   | УПЪТВАНЕ ЗА УПОТРЕБА       | 51 |
| <b>RUS</b>  | ИНСТРУКЦИЯ ПО ЭКСПЛУАТАЦИИ | 54 |
| <b>PL</b>   | INSTRUKCJA OBSŁUGI         | 57 |
| <b>ARAB</b> | دليل الاستخدام             | 60 |
| <b>EST</b>  | KASUTUSJUHEND              | 63 |
| <b>LV</b>   | LIETOŠANAS INSTRUKCIJA     | 66 |
| <b>LT</b>   | NAUDOJIMO INSTRUKCIJA      | 69 |
| <b>RO</b>   | INSTRUCȚIUNI DE UTILIZARE  | 72 |
| <b>SLO</b>  | NAVODILA ZA UPORABO        | 75 |
| <b>SR</b>   | UPUTSTVO ZA UPOTREBU       | 78 |
| <b>HR</b>   | UPUTE ZA UPORABU           | 81 |

**D** Interpretation der Messwerte • **GB** Interpretation of measured values • **NL** Interpretatie van de meetwaarden • **I** Interpretazione dei valori • **F** Interprétation des valeurs de mesure • **E** Interpretación de los valores medidos • **P** Interpretação dos valores medidos • **DK** Tolkning af målingerne • **S** Interpretation av mätvärden • **NO** Tolkning av måleverdiene • **FI** Mittausarvojen tulkinta • **CZ** Interpretácia výsledkov pomiaru • **SK** Interpretácia meraných hodnôt • **BG** Интерпретиране на измерените стойности • **RUS** Интерпретация измерений • **PL** Interpretacja namierzonych wartości • **ARAB** شرح قيم القياسات • **EST** Vaartuse tähendus • **LV** Mērījumu skaidrojums • **LT** Matavimų vertčių paaiškinimas • **RO** Interpretarea valorilor măsurate • **SLO** Razlaga merilnih vrednosti • **SR** Interpretacija mernih vrednosti • **HR** Tumačenje mjernih vrijednosti

| Display    | 1       | 2           | 3           | 4           | 5           | 6           | 7           | 8           | 9        |
|------------|---------|-------------|-------------|-------------|-------------|-------------|-------------|-------------|----------|
|            | —       |             |             | ☺           |             |             | +           |             |          |
| <b>BMI</b> | < 17,4  | 17,5 - 17,9 | 18,0 - 18,4 | 18,5 - 20,4 | 20,5 - 22,9 | 23,0 - 24,9 | 25,0 - 27,4 | 27,5 - 29,9 | ≥ 30     |
| ☹          | < -12 % | - 12 %      | - 7 %       | -3 %        | 0 %         | + 3 %       | + 7 %       | + 12 %      | > + 12 % |
| 💧          | < - 6 % | - 6 %       | - 4 %       | -2 %        | 0 %         | + 2 %       | + 4 %       | + 6 %       | > + 6 %  |
| 🧑          | < - 6 % | - 6 %       | - 4 %       | -2 %        | 0 %         | + 2 %       | + 4 %       | + 6 %       | > + 6 %  |

| Alter/Age | ♀    |         |         |      |        |      | ♂    |         |         |      |        |      |
|-----------|------|---------|---------|------|--------|------|------|---------|---------|------|--------|------|
|           | ☹    |         |         | 💧    |        |      | ☹    |         |         | 💧    |        |      |
|           | —    | ☺       | +       | ++   | ☺      | ☺    | —    | ☺       | +       | ++   | ☺      | ☺    |
| 17-18     | < 15 | ≥ 15-25 | ≥ 25-33 | > 33 | > 52,5 | > 28 | < 8  | ≥ 8-18  | ≥ 18-24 | > 24 | > 58,5 | > 38 |
| 18-30     | < 20 | ≥ 20-29 | ≥ 29-36 | > 36 | > 49,5 | > 28 | < 8  | ≥ 8-18  | ≥ 18-24 | > 24 | > 57   | > 38 |
| 30-40     | < 22 | ≥ 22-31 | ≥ 31-38 | > 38 | > 47,5 | > 28 | < 11 | ≥ 11-20 | ≥ 20-26 | > 26 | > 56   | > 38 |
| 40-50     | < 24 | ≥ 24-33 | ≥ 33-40 | > 40 | > 46,5 | > 28 | < 13 | ≥ 13-22 | ≥ 22-28 | > 28 | > 55   | > 38 |
| 50-60     | < 26 | ≥ 26-35 | ≥ 35-42 | > 42 | > 46   | > 28 | < 15 | ≥ 15-24 | ≥ 24-30 | > 30 | > 54   | > 38 |
| 60+       | < 28 | ≥ 28-37 | ≥ 37-47 | > 47 | > 45   | > 28 | < 17 | ≥ 17-26 | ≥ 26-34 | > 34 | > 53   | > 38 |

**D** Alle Angaben in % • **GB** All figures in % • **NL** Alle waarden in % • **I** Tutte le indicazioni sono in % • **F** Toutes les données en % • **E** Todos los datos en % • **P** Todas as informações em % • **DK** Alle tal i % • **S** Alla uppgifter i % • **NO** Alle tall i % • **FI** Kaikki tiedot prosentteina • **CZ** Wszystkie dane w % • **SK** Všetky údaje sú v % • **BG** Всички данни в % • **RUS** Все показатели в % • **ARAB** % كافة البيانات بالنسبة المئوية • **EST** Kõik väärtused %-na • **LV** Visi lielumi norādīti % • **LT** Visi skaičiai pateikiami procentais (%) • **RO** Toate datele în % • **SLO** Vsi podatki v % • **SR** Svi navodi su u % • **HR** Svi podaci u %

## BMI

**D** Body-Mass-Index • **GB** Body Mass Index • **NL** Body-Mass-Index • **I** Indice di massa corporea (BMI) • **F** Indice de masse corporelle • **E** Índice de masa corporal • **P** Índice de massa corporal • **DK** Body Mass Index • **S** Body-Mass-Index • **NO** Kroppsmasseindeks • **FI** Painoindeksi • **CZ** Wskaźnik masy ciała • **SK** Index telesnej hmotnosti • **BG** Body Mass Index • **RUS** Body-Mass-Index • **PL** Index telesnej hmotności • **ARAB** مؤشر كتلة الجسم • **EST** Kehamassiindeks • **LV** Ķermeņa masas indekss • **LT** kūno masės indeksas (KMI) • **RO** Indexul de masă corporală • **SLO** Indeks telesne mase • **SR** Body Mass Index • **HR** Indeks tjelesne mase

**D** Körperfettanteil • **GB** body fat content • **NL** lichaamsvetaandeel • **I** percentuale di grasso • **F** taux de graisse • **E** porcentaje de la grasa corporal • **P** percentagem de gordura corporal • **DK** fedtprocent • **S** andelen kroppsfett • **NO** Fettprosent • **FI** Kehon rasvaprosentti • **CZ** podíl tělesného tuku • **SK** zawartość tłuszczu w ciele • **BG** Процент на телесната мазнина • **RUS** процент жира в организме • **PL** test zsiřtartalmát • **ARAB** نسبة الدهون في الجسم • **EST** Rasvaprosent • **LV** Ķermeņa tauku procents • **LT** kūno riebalų procentinė dalis • **RO** Coeficient de grăsime corporală • **SLO** Delež telesne maščobe • **SR** Udeo telesne masnoće • **HR** Udio tjelesne masti


**D** Körperwasseranteil • **GB** body water content • **NL** lichaamswateraandeel • **I** percentuale acqua corporei • **F** taux d'eau dans le corps • **E** porcentaje del agua corporal • **P** percentagem de água corporal • **DK** vandprocent • **S** andelen kroppsvatten • **NO** Vannprosent • **FI** Kehon vesipitoisuus • **CZ** podíl vody v těle • **SK** zawartość wody w ciele • **BG** Процент на телесната вода • **RUS** процентное содержание воды в организме • **PL** test wżtartalmát • **ARAB** نسبة الماء في الجسم • **EST** Veeprotsent • **LV** Ķermeņa ūdens procents • **LT** kūno vandens procentinis • **RO** Coeficient de apă în corp • **SLO** Delež vode v telesu • **SR** Udeo telesne vode • **HR** Udio vode u tijelu


**D** Muskelanteil • **GB** muscle proportion • **F** pourcentage de muscles • **I** massa muscolare • **NL** Spierenaandeel • **E** porción de los músculos • **P** percentagem de músculo • **DK** muskler andelen • **S** muskelprosent • **NO** Muskelprosent • **FI** Lihasprosentti • **CZ** podíl svalů • **SK** masy mišičnovej • **BG** Процент на мускулите • **RUS** процент мышц • **PL** izom arány • **ARAB** نسبة العضلات • **EST** Lihasmassi osakaal • **LV** Muskulų procents • **LT** raumenų procentinė dalis • **RO** Coeficient de masă musculară • **SLO** Delež mišične mase • **SR** Udeo mišića • **HR** Udio mišića

# A

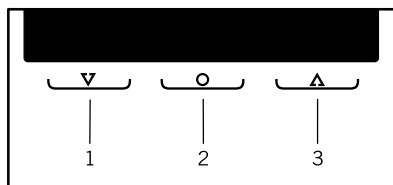
# B
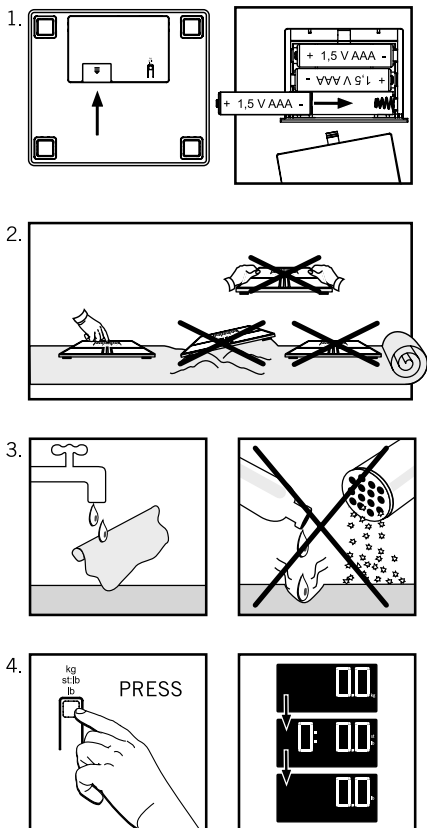
# E


# C

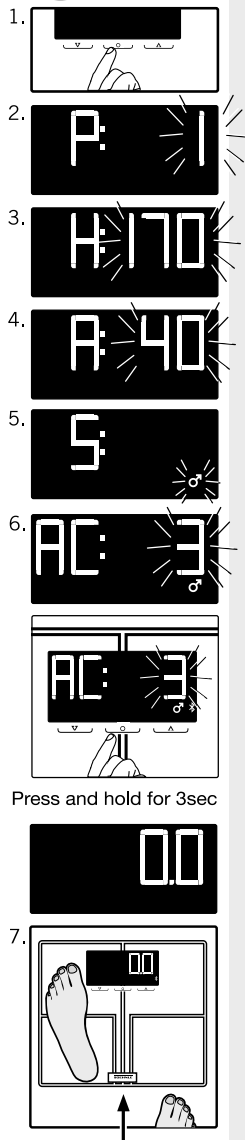Press and hold for 3sec

# D

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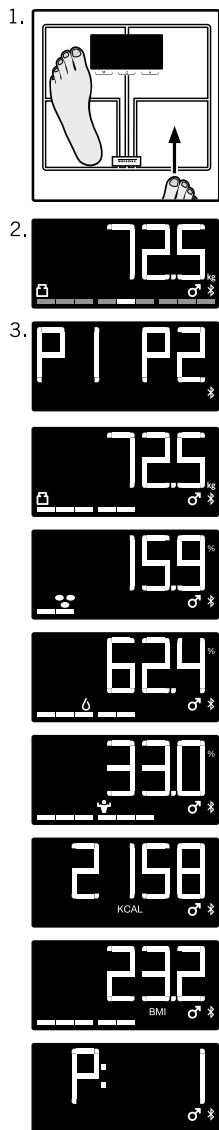**GB** This device may be used by children over the age of 8 and persons with reduced physical, sensory, or mental abilities or with a lack of experience and knowledge only if such are supervised or have been instructed in the safe use of the device as well as the resulting risks. Cleaning and maintenance shall not be made by children without supervision, Children shall not play with the device.

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**Deutschland** Tel.: (08 00) 5 34 34 34

**International** Tel.: +49 26 04 97 70

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Apple, iPhone, iPad, iOS, Apple Health and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries.

Android, Google Fit, Google Play and the Google Play logo are trademarks of Google Inc.

The **Bluetooth®** word mark and logos are registered trademarks owned by Bluetooth SIG, Inc.

## A. Operating elements

1. Minus (▼)
2. Confirm (○)
3. Plus (▲)

Sensor touch:

Only touch the keys softly!

## B. Preparation

1. Insert batteries (3 x 1.5 V AAA).

### Note:

**Immediately** (as long as "0.0" is still being displayed) place the body analysis scale on a level surface and wait until the scale powers down automatically.

Then begin with data entry. Otherwise, it might be possible that the weight of the scale would erroneously be added to the initial weighting process. The same might occur during regular operation if the scale is held in the hand prior to the weighting process.

2. For all measurements, place scale level on a hard surface (not on carpet).
3. Cleaning and care: Clean only using a slightly damp cloth. Do not use solvents or abrasives. Never submerge scale in water.
4. A country-specific option to switch from kg/cm to st/in or lb/in is possible by pushing the key on the bottom of the scale to switch units

**Attention!** Slipping hazard in case of wet surface.

## C. Data input

1. Press ○ button.
2. Use ▼/▲ buttons to select memory slot and confirm with ○. → P
3. Use the ▼/▲ buttons to set your height and confirm with ○. → H
4. In the next step, set your age and confirm with ○. → A
5. Set gender and confirm with ○. → S
6. To finish, set the activity level and confirm with ○. (Hold down key for 3 sec.). → AC

### Athlete mode:

By setting the personal activity level, the respective calorie requirement can be precisely determined.

| AC | Explanation                                                       | Hrs per week                            |
|----|-------------------------------------------------------------------|-----------------------------------------|
| 1  | Hardly any physical activity                                      |                                         |
| 2  | Slight, active, sitting and standing activity                     |                                         |
| 3  | House and garden work, mainly standing, sitting only occasionally |                                         |
| 4  | Athletes and people with plenty of movement                       | Intensive training 5-10 hrs/week        |
| 5  | Heavy manual labourer, highly active athletes                     | Intensive training at least 10 hrs/week |

7. Now **immediately** place the body analysis scale on the floor, wait until the screen displays "0.0" and then step on the scale with bare feet.

The initial measurement is indispensable and serves to store the data for later automatic recognition of persons. If no person steps on the scale, data input must be repeated. The first body analysis is performed after a person steps on the scale.

This scale is equipped with a Time-Out mode. It will automatically switch off after approx. 40 seconds if no keys are pressed.

### Alternative:

You can also perform this step conveniently via the Soehnle MyScale App, in which case the data is transferred to the scale. The first measurement is necessary, however. Once you have entered the data on the scale for the recognition of persons as described under C, the Soehnle MyScale App will automatically adopt this.

## D. Body analysis

1. Place scale on a solid, level floor and wait until scale switches off automatically.  
Step on the scale with bare feet. Stand still.
2. Body weight and the other calculated data are displayed.
3. If the calculated values cannot be clearly allocated to one user, the possible memory slots will be displayed.

If the calculated weight can be allocated to more than 2 users, other possible memory slots will be displayed after 3 seconds.

Use the ▼/▲ keys (▼= left, ▲= right) to confirm the memory slot displayed there.

The scale measures the individual body values by means of bioimpedance analysis (body fat percentage, body water, muscle mass), calculates the basic metabolic rate per day (kcal) and the BMI (Body Mass Index). The 9-segment status bar shows the respective deviation from the normal value (see overleaf).

The scale switches off automatically.

**Note:** If needed, the user can step off the scale as soon as the first value appears on the display.

### Interpretation of measured values

In order to be better able to assess your measured values, at the front of this operating manual you will find an overview table, from which you can take the commonly accepted optimal values.

The optimal range is marked in the table with a ☺. Values outside this range are either too low (–) or too high (+/++).

Here is an example:

A female user is 35 years old and is shown a body fat content of 28.2 % from the scale. For the age group of 30-40 year old women, this is a value in the optimal range (22-31 % → ☺).

### Classification of weight according to BMI (source: World Health Organisation, 2016)

| BMI (kg/m <sup>2</sup> ) | Classification    |
|--------------------------|-------------------|
| < 18,5                   | Underweight       |
| 18,5 - 24,9              | Normal weight     |
| 25,0 - 29,9              | Overweight        |
| 30,0 - 34,9              | Obesity level I   |
| 35,0 - 39,9              | Obesity level II  |
| > 40,0                   | Obesity level III |

## E. Notifications

1. Replace battery (Lo = low battery)
2. Overloading of scale (max. 180 kg)
3. The measured value is above (H) or below (L) the threshold value specified below.

|   | L  H | L  H | L  H |
|---|-----------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------|
| ♀ | < 10 %   > 60 %                                                                         | < 5 %   > 90 %                                                                          | < 5 %   > 90 %                                                                          |
| ♂ | < 4 %   > 60 %                                                                          | < 5 %   > 90 %                                                                          | < 5 %   > 90 %                                                                          |

## F. Installation of the Soehnle MyScale App

Install the Soehnle MyScale App on your smartphone or tablet as follows:

### Android™ smartphones / tablets

If you use an Android device, download and install the Soehnle MyScale App from the Google Play Store. To do so, scan the QR code or open the app with the name "Play Store" and look for "Soehnle MyScale App". Then load and install the free app.

### iPhone® / iPad® (iOS®)

If you use an Apple iPhone or iPad, download and install the Soehnle MyScale App from the Apple App Store. To do so, scan the QR code or open the app with the name "App Store" on your iPhone and look for "Soehnle MyScale App". Then load and install the free app.

### Compatibility

The Soehnle MyScale App is compatible with Google Fit® and Apple Health®.

## G. Note on connecting the scale to a smartphone/tablet

The scale can only be connected to a smartphone or tablet using the Soehnle MyScale app. Please make sure that **Bluetooth** is switched on on your device. It is not possible to pair the scale via the Android™ or iOS® operating system itself. The connection can only be established via the Soehnle MyScale app.

## H. Technical data

**Bluetooth®:** 4.1

Frequency band: 2.4 GHz

Maximum transmission power: 10 mW

Max. 180 kg / 100 g

Max. 31 st / 0.2 lb

Max. 396 lb / 0.2 lb

Age: 17-99 years

Batteries required: 3 x 1.5 V AAA

### Guarantee

Leifheit AG offers a 3 year guarantee for this quality product as of the purchase date (or in case of orders, as of receipt of the goods). Guarantee claims must be asserted within the guarantee period immediately after occurrence of the defect. The guarantee includes product quality.

The following are excluded from the guarantee:

- (1) Wear caused by use or other natural wear;
- (2) Damage caused by improper use or handling (e.g. impact, shock, drop),
- (3) Damage caused by non-compliance with these operating instructions;
- (4) Battery, or rechargeable battery

In case of a guarantee claim, Leifheit, at its own discretion, either repairs the defective components or replaces the product. If a repair cannot be performed and an identical product is no longer available in the product range for the purpose of replacement, you will receive a replacement product of the closest equivalent value. Reimbursement of the purchase price is not possible in case of a guarantee claim. This guarantee furthermore does not grant any claims for damages.

To assert the guarantee, please contact the seller where you have acquired the product and present the defective product and the purchase receipt (copy). This guarantee is valid worldwide.

Your statutory rights, in particular warranty rights, continue to apply and are not limited by this guarantee.

### CE conformity

This device complies with the applicable EU Directive ([www.soehnle.com](http://www.soehnle.com)).

### Battery disposal EU Directive

Batteries are not a part of your regular household waste. You must return batteries to your municipality's public collection or wherever batteries of the respective type are being sold.

Batteries that contain hazardous substances feature this sign:

Pb = Battery contains lead

Cd = Battery contains cadmium

Hg = Battery contains mercury

### Disposal of electric and electronic devices EU Directive

This product is not to be treated as regular household waste but must be returned to a collection point for recycling electric and electronic devices. Further information is available from your municipality, your municipality's waste disposal services, or the retailer where you purchased your product.

### Consumer service

The following contact partners are happy to help you in case of questions and comments:

**Germany** Tel.: (08 00) 5 34 34 34

**International** Tel.: +49 26 04 97 70

**Monday to Thursday** 8.00 am to 5.00 pm

**Friday** 8.00 am to 4.00 pm

### Legal information

Apple, iPhone, iPad, iOS, Apple Health and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries.

Android, Google Fit, Google Play and the Google Play logo are trademarks of Google Inc.

The **Bluetooth®** word mark and logos are registered trademarks owned by Bluetooth SIG, Inc.

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The **Bluetooth**\* word mark and logos are registered trademarks owned by Bluetooth SIG, Inc.
Android, Google Fit, Google Play i logotip trgovine Google Play su žigovi tvrtke Google Inc.
Verbalni žig **Bluetooth®** i logotipovi registrirani su žigovi u vlasništvu tvrtke Bluetooth SIG, Inc.