SOEHNLE

Shape Sense Control 200

D	BEDIENUNGSANLEITUNG	2
GB	USER MANUAL	4
F	MODE D'EMPLOI	6
I	ISTRUZIONI PER L'USO	8
NL	GEbruIKSAANWIJZING	10
E	INSTRUCCIONES DE MANEJO	12
P	MANUAL DE INSTRUÇÕES	14
DK	BRUGSANVISNING	16
S	BRUKSANVISNING	18
CZ	INSTRUKCJA OBSŁUGI	20
SK	NÁVOD NA OBSLUHU	22
PL	NÁVOD K POUŽITÍ	24
RUS	ИНСТРУКЦИ ПО ЭКСПЛУАТАЦИИ	26
TR	KULLANMA KILAVUZU	28
LT	NAUDOJIMO INSTRUKCIJA	30
BG	УПЪТВАНЕ ЗА УПОТРЕБА	32
ARAB	دليل الاستخدام	34

- Nicht geeignet für Personen mit elektronischen Implantaten (z.B. Herzschrittmacher).
- It is not suitable for persons with electronic implants (pacemaker, etc.).
- Les personnes qui portent des implants médicaux (p. ex. un stimulateur cardiaque etc.) ne doivent pas utiliser cet appareil.

- Non è adatto per persone con impianti elettronici (pacemaker, ecc.).
- Niet geschikt voor personen met elektronische implantaten (pacemaker, enz.).
- No es apropiado para personas con implantes electrónicos (marcapasos, etc.).
- Não é adequado para pessoas com implantes eletrônicos (pacemakers, etc.).
- Ej lämplig för personer med elektroniska implantat (pace-maker etc).
- Ikke egnet til folk med elektroniske implantater (pacemakere ect).
- Nie nadaje się dla osób, którym wszczepiono elektroniczne implanty (np. rozruszniki serca).
- Pístroj nie je vhodný pre osoby s elektronickými implantátmi (napr. kardiostimulátorom).
- Není vhodný pro osoby s elektronickými implantáty (kardiostimulátory, atd.).
- Никогда не используйте данный прибор вместе с электронными медицинскими приборами, например, кардиостимуляторами.
- Elektronik parçalar (kalp pili, vs.) taşıyan insanlar için uygun değildir.
- Netinka elektroninius implantus (pvz. širdies stimulatorius) turintiems žmonėms.
- Уредът не е подходящ за хора с електронни импланти (напр. сърдечни пейсмейкери).
- غير مناسب للأشخاص الذين يستخدمون أجهزة إلكترونية مزروعة (منظم ضربات القلب).

Alter/Age	♀						♂					
	-		+	++			-		+	++		
17-18	< 15	≥ 15-25	≥ 25-33	> 33	> 52,5	> 28	< 8	≥ 8-18	≥ 18-24	> 24	> 58,5	> 38
18-30	< 20	≥ 20-29	≥ 29-36	> 36	> 49,5	> 28	< 8	≥ 8-18	≥ 18-24	> 24	> 57	> 38
30-40	< 22	≥ 22-31	≥ 31-38	> 38	> 47,5	> 28	< 11	≥ 11-20	≥ 20-26	> 26	> 56	> 38
40-50	< 24	≥ 24-33	≥ 33-40	> 40	> 46,5	> 28	< 13	≥ 13-22	≥ 22-28	> 28	> 55	> 38
50-60	< 26	≥ 26-35	≥ 35-42	> 42	> 46	> 28	< 15	≥ 15-24	≥ 24-30	> 30	> 54	> 38
60+	< 28	≥ 28-37	≥ 37-47	> 47	> 45	> 28	< 17	≥ 17-26	≥ 26-34	> 34	> 53	> 38

D Alle Angaben in % • GB All figures in % • F Toutes les données en % • I Tutte le indicazioni sono in % • NL Alle waarden in % • E Todos los datos en %
 • P Todas as informações em % • DK Alle tal i % • S Alla uppgifter i % • CZ Wszystkie dane w % • SK Všetky údaje sú v % • PL Všechny údaje jsou v % • RUС
 Все показатели в % • TR Tüm rakamlar % cinsindendir • LT Visi skaičiai pateikiami procentais (%) • BG Всички данни в % • ARAB كافة البيانات بالنسبة المئوية %

D Körperfettanteil • GB body fat content • F taux de graisse • I percentuale di grasso • NL lichaamsvetaandeel
 • E porcentaje de la grasa corporal • P percentagem de gordura corporal • DK fedtprocent • S andelen kroppsfett • CZ test
 zsrírtartalmát • SK zawartość tłuszczu w ciele • PL podіл тілесного тuku • RUС процент жира в организме • TR vücudtaki yağ oranını
 • LT Kūno riebalų procentinė dalis • BG Процент на телесната мазнина • ARAB نسبة الدهون في الجسم



D Körperwasseranteil • GB body water content • F taux d'eau dans le corps • GB percentuale acqua corporei • NL lichaams-
 wateraandeel • E porcentaje del agua corporal • P percentagem de água corporal • DK vandprocent • S andelen kroppsvatten
 • CZ test vizírtartalmát • SK zawartość wody w ciele • PL podіл vody v těle • RUС процентное содержание воды в организме
 • TR ücudtaki su oranını • LT Kūno vandens procentinis • BG Процент на телесната вода • ARAB نسبة الماء في الجسم







D Muskelanteil • GB muscle proportion • F pourcentage de muscles • I massa muscolare • NL Spierenaandeel
 • E porción de los músculos • P percentagem de músculo • DK muskler andelen • S muskelen • CZ izom arány,
 masy mięśniowej • SK podílú svalů • RUС процент мышц • TR vücudtaki kas oranının • LT raumenų procentinė dalis
 RUС Процент на мускулите • ARAB نسبة العضلات





A




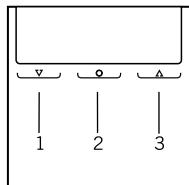
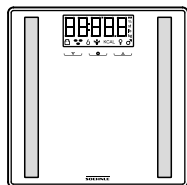
B

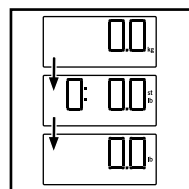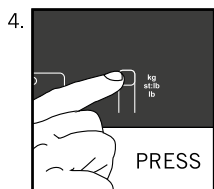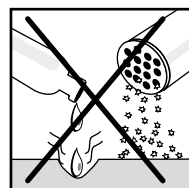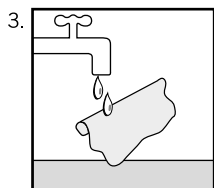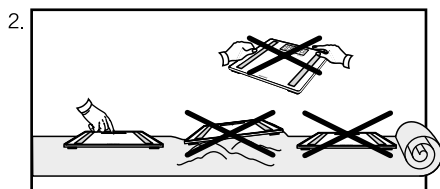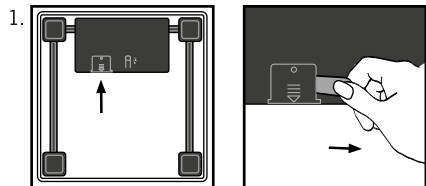C

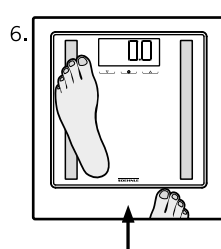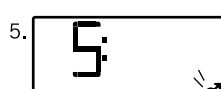
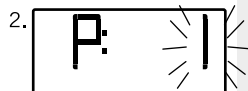
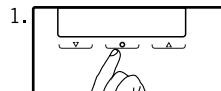
D

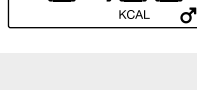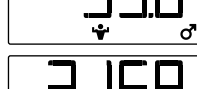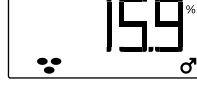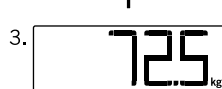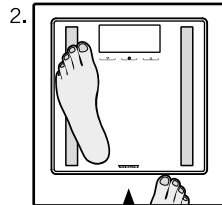E

New battery 1 x 3 V CR2430

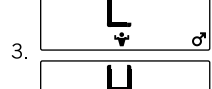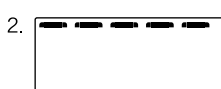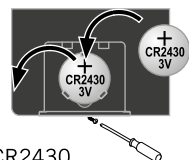A. Bedienelemente

1. Minus (▼)
2. Bestätigen (○)
3. Plus (▲)

Sensor-Touch:

Tasten nur sanft berühren!

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1. Batteriestrip entfernen.

Hinweis:

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Touch keys only slightly!

B. Preparation

- 1. Remove battery strip.
Note: Immediately (as long as „0.0“ is still being displayed) place the body analysis scale on a level surface and wait until the scale powers down automatically.
Then begin data input. Otherwise, it might be possible that the weight of the scale would erroneously be added to the initial weighting process.
The same might occur during regular operation if the scale is held in the hand prior to the weighting process.
- 2. For all measurements, place scale level on a hard surface (not carpet).
- 3. Cleaning and Care: clean only using slightly damp cloth. Do not use solvents or abrasives. Never submerge scale in water.
- 4. A country-specific option to switch from kg/cm to lb/in or st/in is possible by pushing the key on the bottom of the scale to switch units.

Attention! Slipping hazard in case of wet surface.

C. Data input

- 1. Push ○ key.
- 2. Use ▲/▼ keys to select memory slot and confirm (○).
- 3. Then use ▲/▼ keys enter body height and confirm.
Note: For unit conversion, see section B.4.
- 4. In the next step enter the age and confirm.
- 5. Finally select sex and confirm (Press ○ key 3 sec.).
- 6. Now immediately place the body analysis scale on the floor, wait until the screen displays „0.0“ and then step on the scale with bare feet.
The initial measurement is indispensable and

serves to store the data for later automatic recognition of persons. If no person steps on the scale, data input must be repeated.
The first body analysis is performed after a person steps on the scale.
This scale is equipped with a Time-Out mode. It will automatically switch off after approx. 40 seconds if no keys are pressed.

Resetting the scale:

With the scale switched off, press and hold the unit button (B4) on the back of the scale for approximately 5 seconds until **dEL** is displayed and confirm with the ○ -button – all data will be deleted and the scale reset to factory settings.

D. Body analysis

- 1. Place scale on hard, level surface and wait until the scale switches off automatically.
- 2. Step barefoot onto scale. Please stand still.
- 3. The body weight will be displayed.
The scale determines the body values (body fat, body water, muscle mass, kcal).
Note: Body values are calculated based on the determined weight and the personal data.
The scale shuts off automatically.

E. Notifications

- 1. Replace batteries
- 2. Overload: max. 180 kg
- 3. The measured value is higher (H) or lower (L) than the limit.

	L ☼ H		L 💧 H		L 🧑 H	
♀	< 10 %	> 60 %	< 5 %	> 90 %	< 5 %	> 90 %
♂	< 4 %	> 60 %	< 5 %	> 90 %	< 5 %	> 90 %

F. Technical data

Max. 180 kg / 100 g
Max. 396 lb / 0.2 lb
Max. 28 st / 0.2 lb
Age: 17-99 years
Battery requirement: 1 x 3V CR 2430

Guarantee

Leifheit AG grants you 3 years guarantee as of the date of purchase (or for order upon receipt of the goods) for the quality product at hand. The guarantee claims must be asserted immediately upon occurrence of the defect within the guarantee period. The guarantee is limited to the quality of the products.

The following are excluded from the guarantee:

- (1) wear caused by proper use or natural causes.
- (2) damage caused by improper use or handling (e.g. impact, bump, fall).
- (3) damage caused by non-compliance with stated operating instructions.
- (4) the rechargeable battery.

In case of a guarantee claim, Leifheit will at its own discretion provide either repair of defective parts or replace the product. If a repair is not possible and if an identical product is no longer available within the product range for replacement, you will receive a replacement product that is, insofar as possible, of equal value. In case of a guarantee claim, the purchase price cannot be reimbursed. This guarantee also does not grant compensation rights. This guarantee furthermore does not grant any claims for damages.

To assert the guarantee, please contact the seller where you purchased the product and present the defective product, the purchase receipt (copy) and the guarantee slip.

The guarantee is valid globally.

Your statutory rights, in particular warranty rights, continue to be valid and are not limited by this warranty.

CE conformity

This device complies with the applicable EU Directive (www.soehnle.com).

Battery disposal EU Directive

Batteries are not a part of your regular household waste. You must return batteries to your municipality's public collection or wherever batteries of the respective type are being sold.

Batteries containing hazardous materials are marked with these symbols:

Pb = Battery contains lead

Cd = Battery contains cadmium

Hg = Battery contains mercury

Disposal of electric and electronic devices EU Directive

This product is not to be treated as regular household waste but must be returned to a collection point for recycling electric and electronic devices. Further information is available

from your municipality, your municipality's waste disposal services, or the retailer where you purchased your product.

Consumer service

If you have questions or comments, the following points of contact are available:

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